

Ranitidine(Degraded) by HPLC - Luna 3µm C8(2) 100x4.6mm

Column: Luna® 3 µm C8(2) 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4248-E0

Elution Type: Gradient

Eluent A: 98:2 40mM potassium phosphatate pH 7.1:Acetonitrile

Eluent B: 78:22 40mM potassium phosphatate pH 7.1:Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	10	0	100
	3	15	0	100
	4	16	100	0
	5	20	100	0

Flow Rate: 1.5 mL/min

Col. Temp.: 35 °C


Detection: UV-Vis Abs.-Diode Array (PDA) @ 230 amu (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

19506

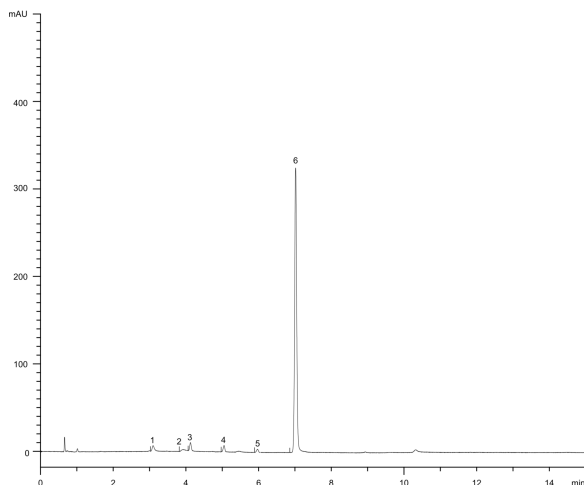

ANALYTES:

1 Ranitidine amino alcohol hemifumarate

Retention Time: 3.096 min

2 Ranitidine diamine hemifumarate

Retention Time: 3.926 min

3 Ranitidine S-oxide

Retention Time: 4.12 min

4 Ranitidine N-oxide

Retention Time: 5.043 min

5 Ranitidine complex nitroacetamide

Retention Time: 5.969 min

6 Ranitidine

Retention Time: 7.018 min



Products used in this application:

