

ABN16 (acids, bases, neutrals) Mix using Kinetex 2.6µm C8, 50x2.1mm

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4497-AN

Elution Type: Gradient

Eluent A: 0.1% formic acid/water

Eluent B: 0.1% formic acid/ACN


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.2	95	5
	3	4.2	5	95
	4	4.21	95	5
	5	5.5	95	5

Flow Rate: 0.8 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19539

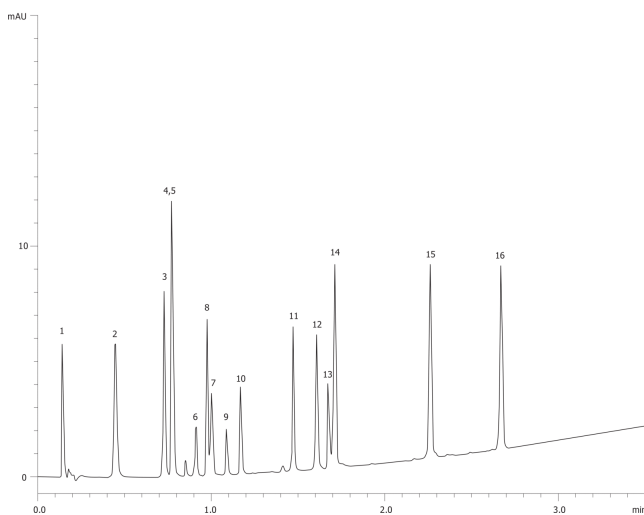

Products used in this application:



ABN16 (acids, bases, neutrals) Mix using Kinetex 2.6µm C8, 50x2.1mm

ANALYTES:

- 1** Pyridine
- 2** Acetaminophen
- 3** Sulfathiazole
- 4** Pindolol
- 5** Quinidine
- 6** Benzyl alcohol
- 7** Acebutolol
- 8** Phenol
- 9** Chlorpheniramine
- 10** Triprolidine
- 11** Prednisolone
- 12** 3-Methyl-4-nitrobenzoic acid
- 13** Nortriptyline
- 14** 2-Hydroxy-5-methylbenzaldehyde
- 15** Diflunisal
- 16** Hexanophenone

