

Prodigy Test Mix using Kinetex 2.6µm C8 50x2.1mm ID

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4497-AN

Elution Type: Isocratic

Eluent A: Water

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19555



Products used in this application:

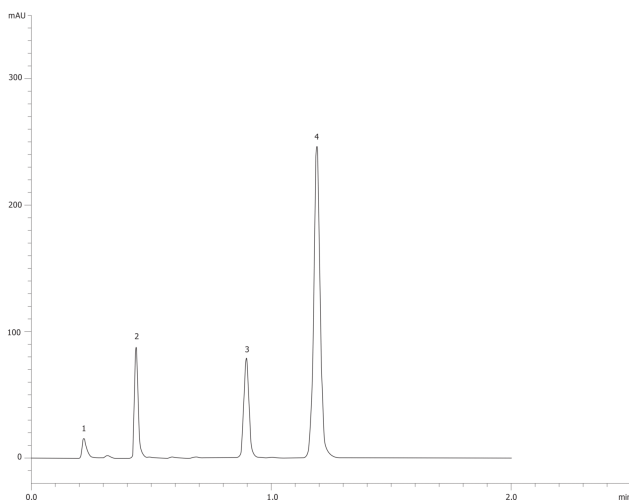

ANALYTES:

1 Uracil

Retention Time: 0.216 min

2 Acetophenone

Retention Time: 0.402 min

3 Toluene

Retention Time: 0.762 min

4 Naphthalene

Retention Time: 0.959 min

