

USP Loratadine on Kinetex 2.6µm C8 100x4.60mm

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4497-E0

Elution Type: Isocratic

Eluent A: 31.5:31.5:36.8 Methanol:Acetonitrile:Solution A(1.74mg/mL Anhydrous Dibasic Potassium)

Gradient

Step No.	Time (min)	Pct A
----------	------------	-------


Profile:	1	15	100
-----------------	----------	----	-----

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19563

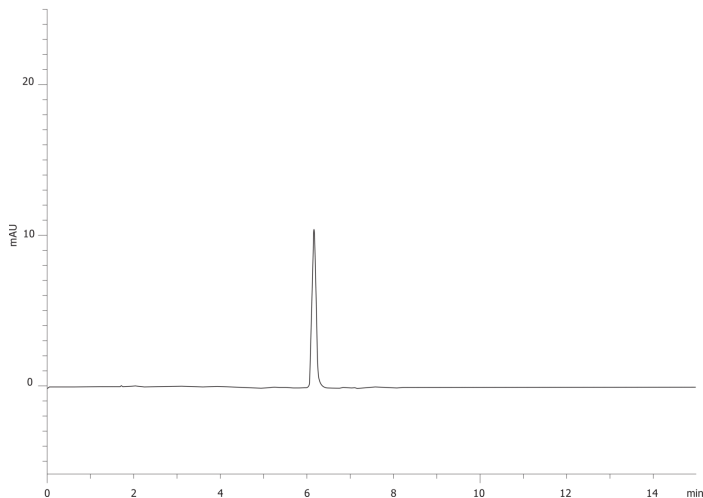

ANALYTES:

1 Loratadine



Products used in this application:

