


USP Loratadine on Kinetex 2.6µm C8 75x4.60mm

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 75 x 4.6 mm, Ea
Dimensions: 75 x 4.6 mm ID
Order No: 00C-4497-E0
Elution Type: Isocratic
Eluent A: 31.5:31.5:36.8 Methanol:Acetonitrile:Solution A(1.74mg/mL Anhydrous Dibasic Potassium)
Gradient

Step No.	Time (min)	Pct A
1	15	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19564

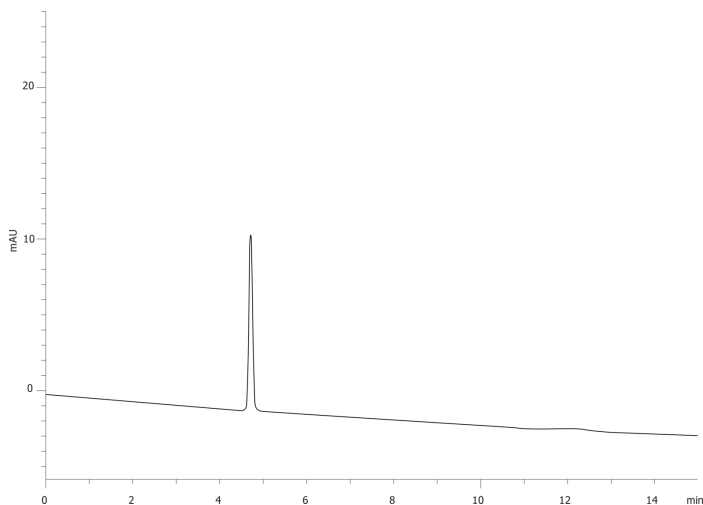

ANALYTES:

1 Loratadine



Products used in this application:

