

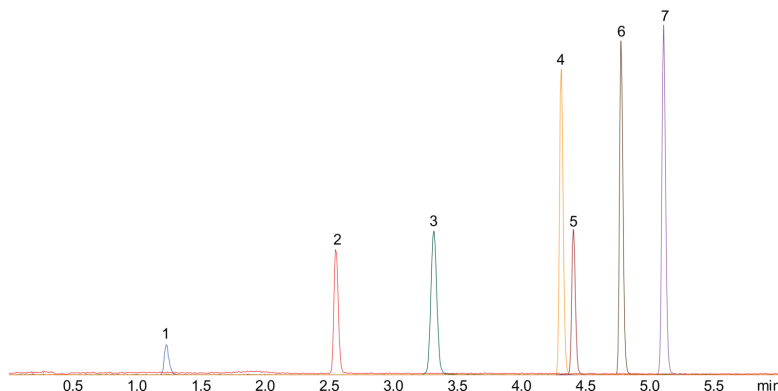
Triazine Pesticides (EPA 536) on Kinetex 2.6µm C8 by LC/MS/MS**Column:** Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4497-AN**Elution Type:** Gradient**Eluent A:** 5 mM Ammonium Acetate**Eluent B:** 100% MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.25	60	40
	3	2	60	40
	4	3	25	75
	5	4	25	75
	6	4.1	95	5

Flow Rate: 0.3 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** <a target="_blank"

19574

ref="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**Products used in this application:****ANALYTES:****1** Atrazine-desethyl-desisopropyl

Retention Time: 1.23 min

2 Atrazine-desisopropyl

Retention Time: 2.55 min

3 Atrazine-desethyl

Retention Time: 3.31 min

4 Cyanazine

Retention Time: 4.3 min

5 Simazine

Retention Time: 4.4 min

6 Atrazine

Retention Time: 4.77 min

7 Propazine

Retention Time: 5.11 min

