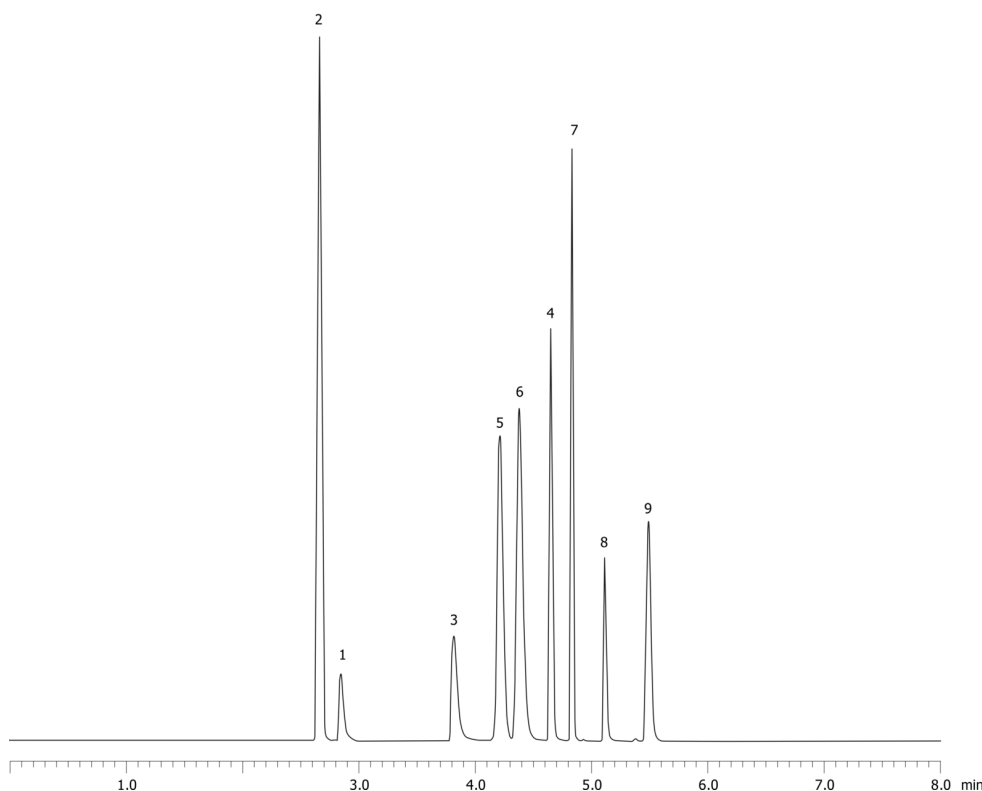


Optimized Separation of Propanal from Blood Alcohols on ZB-BAC2, 30 x 0.53 x 2.00**Column:** Zebron[™] ZB-BAC-2, GC Cap. Column 30 m x 0.53 mm x 2.00 μ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 2 μ m**Order No:** 7HK-G022-32**Oven Profile:** 30 °C Isothermal**Carrier Gas:** Ramp Flow Helium, 4 mL/min**Injection:** Headspace 30:1 1 μ L @ 180°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Sample was 0.100% BAC in water, analyzed by headspace GC
Run as a dual column setup with ZB-BAC1
Ramped flow program was 4 mL/min. for 4.5 min. to 40 mL/min. at 80 mL/min.2 for rest of run.
Liner was 2 mm. i.d. straight**Products used in this application:**

19581

**ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Propanal
- 6 Acetone
- 7 tert-Butanol
- 8 n-Propanol
- 9 2-Butanol

©2005 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com

Phenomenex products are available worldwide.

www.phenomenex.com.cninfo@phenomenex.com