

Phencyclidine from urine using Strata-X-Drug B and Kinetex C18, 50 x 2.1 mm by LCMS

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid


Eluent B: Methanol with 0.1% Formic Acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	1	5	95
	3	2	5	95
	4	2.1	75	25
	5	4	75	25

Flow Rate: 400 µL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ (500 °C)

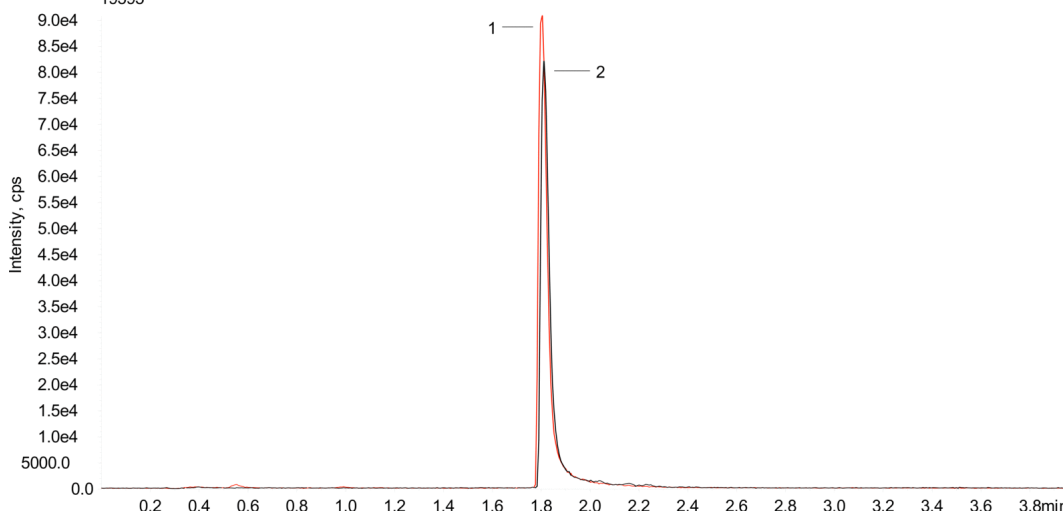

Products used in this application:



ANALYTES:

1 Phencyclidine-d5

Retention Time: 1.82 min

2 Phencyclidine

Retention Time: 1.82 min

