

Kinetex 100 x 2.1 QC with Unmodified 1100

Column: Kinetex® 2.6 µm XB-C18 100 Å, LC Column 100 x 2.1 mm, Ea

Dimensions: 100 x 2.1 mm ID

Order No: 00D-4496-AN

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19614



Products used in this application:

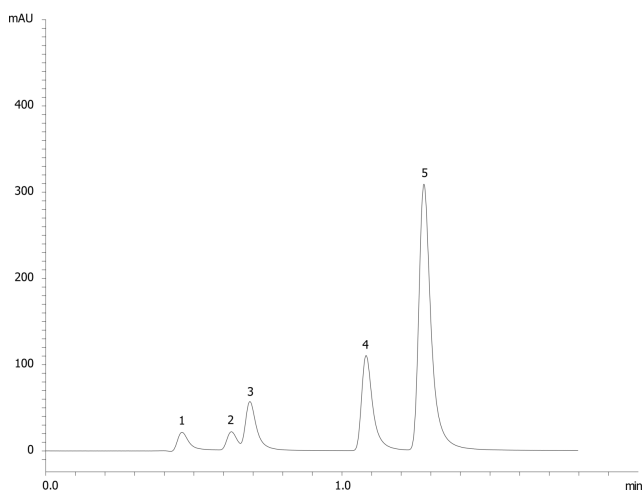

ANALYTES:

- 1** Uracil
Retention Time: 0.46 min
- 2** impurity
Retention Time: 0.59 min
- 3** Acetophenone
Retention Time: 0.689 min
- 4** Toluene
Retention Time: 1.08 min
- 5** Naphthalene
Retention Time: 1.28 min

