

HPLC Application

ID No.: 19684

Ranitidine Tablet on HPLC - Kinetex 2.6µm C18 100x4.6mm - 0.1N HCl Degradation

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Gradient

Eluent A: 98:2 40mM Potassium Phosphate pH 7.1:ACN

Eluent B: 78:22 40mM Potassium Phosphate pH 7.1:ACN

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	100	0
2	10	0	100

Profile:

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 230 nm (ambient)

19684



Products used in this application:

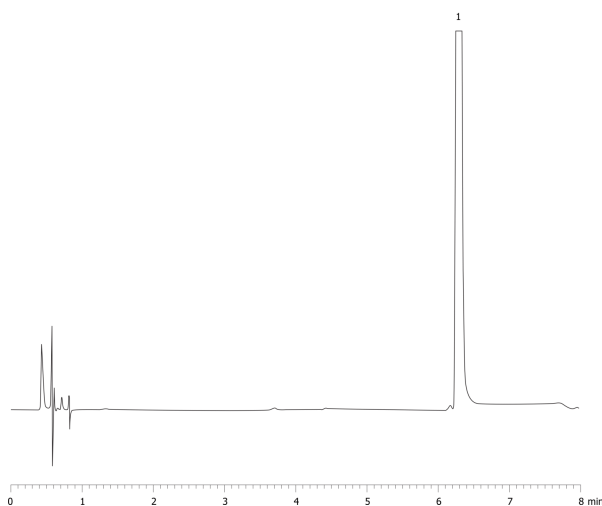

ANALYTES:

1 Ranitidine

Retention Time: 6.19 min

