

## LC/MS/MS Analysis of Digoxin and Digitoxin in plasma (0.25 ng/mL) using Strata-X and Kinetex C8

**Column:** Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4497-AN

**Elution Type:** Gradient

**Eluent A:** 10mM Ammonium acetate

**Eluent B:** 10mM Ammonium acetate in methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50
	2	2.5	0	100
	3	2.51	50	50
	4	5	50	50

**Flow Rate:** 400 µL/min

**Col. Temp.:** ambient

**Detection:** Mass Spectrometer (MS) @ (350 °C)

**Detector Info:** <a target="\_blank"

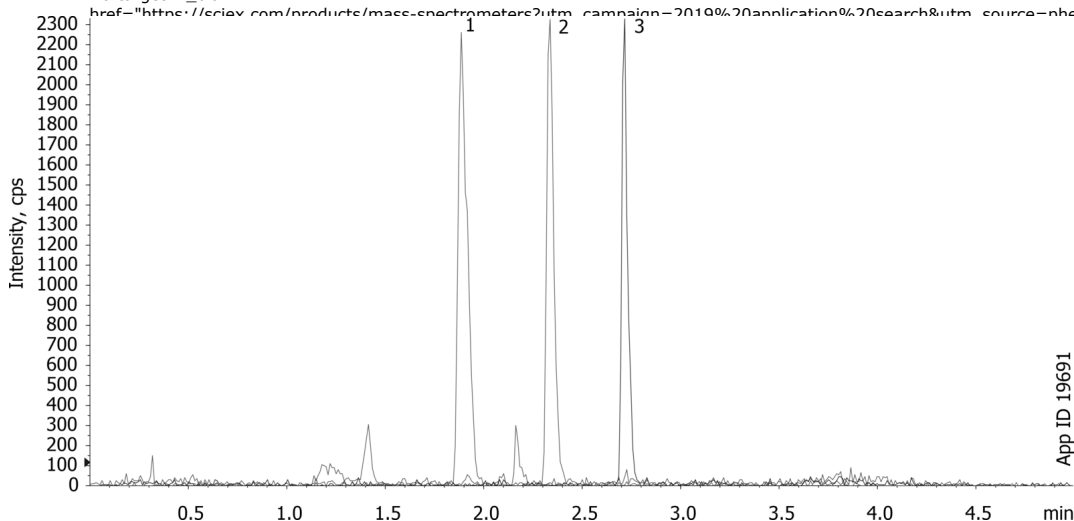

### ANALYTES:

- 1** Digoxin  
Retention Time: 1.9 min
- 2** Oleandrin  
Retention Time: 2.35 min
- 3** Digitoxin  
Retention Time: 2.75 min



Products used in this application:



# Sample Preparation Details

for HPLC Application ID No.: 19691

## LC/MS/MS Analysis of Digoxin and Digitoxin in plasma (0.25 ng/mL) using Strata-X and Kinetex C8

### PRODUCT DESCRIPTION:

Strata™-X 33 µm Polymeric Reversed Phase, 30 mg / 3 mL, Tubes , 50/Pk

Order No.: 8B-S100-TBJ

### SOLID PHASE EXTRACTION (SPE) PROCEDURE:

**Note:** The solvent volumes shown below are for a 30 mg bed mass.

The solvent volumes will need to be adjusted for a smaller or larger bed mass.

#### Condition:

#### Load:

#### Wash:

#### Dry:

10min at full vacuum

#### Elute:

#### Final Prep and Analysis:

Inject: 0 µL on HPLC Mass Spectrometer (MS) @ (350°C)

ANALYTES:	Spiked Conc. (ng/mL)	Log P	pKa	% Rec	%RSC (n=0)
1 Digoxin	0.25			87	
2 Oleandrin	0.25				
3 Digitoxin	0.25			108	

**Note:** This method is designed as a convenient starting point for further investigation and can be tailored to meet your extraction goals.  
Call your local Phenomenex Representative for assistance in method development and optimization techniques.

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at [info@phenomenex.com](mailto:info@phenomenex.com)



Phenomenex products are available worldwide.

[www.phenomenex.com.cn](http://www.phenomenex.com.cn)

[info@phenomenex.com](mailto:info@phenomenex.com)