


Phthalates for China on Zebron ZB-MR-1 20m x 0.18mm x 0.18um**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 20 m x 0.18 mm x 0.18 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 20 meters x 0.18 mm x 0.18 μ m**Order No:** 7FD-G016-08**Oven Profile:** 80C to 320C at 30C/min hold for 5 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Split 10:1 1 μ L @ 260°C**Detection:** Mass Spectrometer (MS) (230°C)**Analyst Note:** Diisononyl phthalate and Diisodecyl phthalate are 400 ppm with all other compounds at 100 ppm in methanol.**Products used in this application:****ANALYTES:**

- 1 Dibutyl phthalate
- 2 Benzyl butyl phthalate
- 3 Di-2-ethyl hexyl phthalate
- 4 Di-n-octyl phthalate
- 5 Diisononyl phthalate
- 6 Diisodecyl phthalate

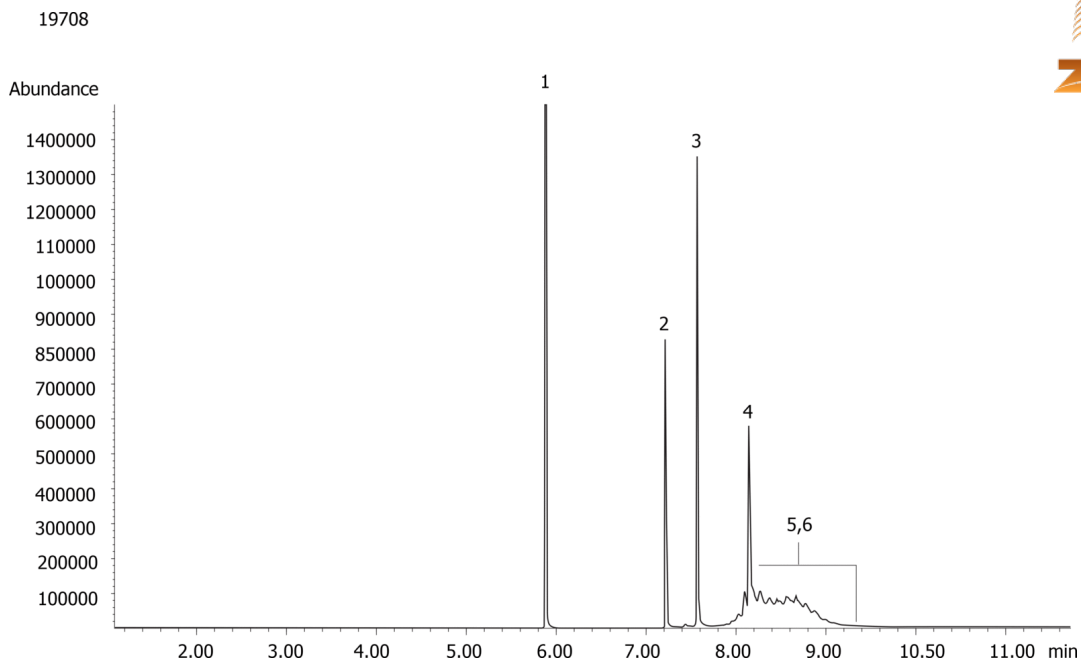