

Benzodiazepines from urine using Strata-X-Drug N on a Kinetex C18 2.6µm, 50x2.1mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid

Eluent B: 0.1% Formic acid in methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35
	2	4	5	95
	3	4.01	65	35
	4	7	65	35

Flow Rate: 300 µL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

Detector Info: <a target="_blank"


href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

Source: Cur-11, CAD-12, IS-4500, Tem-450, NEB-12

Dwell: 25ms

Polarity: Positive

19714

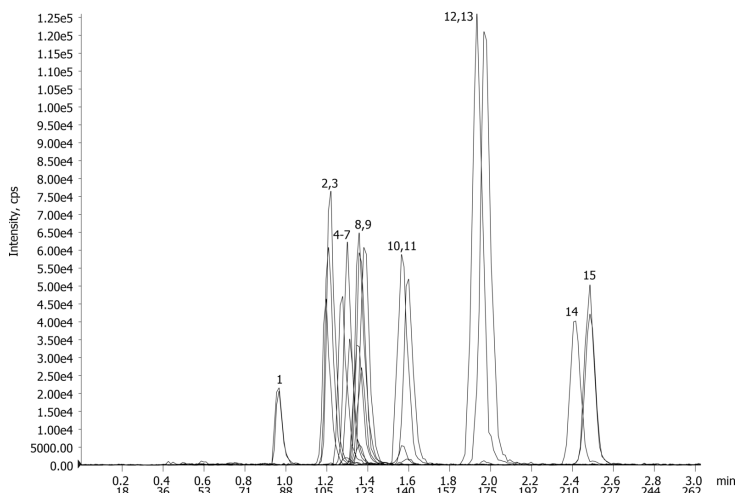

Products used in this application:



Benzodiazepines from urine using Strata-X-Drug N on a Kinetex C18 2.6µm, 50x2.1mm

ANALYTES:

- 1** alpha-Hydroxyalprazolam
- 2** Oxazepam-d5
- 3** Oxazepam
- 4** Alprazolam-d5
- 5** Alprazolam
- 6** Nordiazepam-d5
- 7** Nordiazepam
- 8** Lorazepam-d4
- 9** Lorazepam
- 10** Clonazepam-d4
- 11** Clonazepam
- 12** Temazepam-d5
- 13** Temazepam
- 14** Diazepam-d5
- 15** Diazepam



Sample Preparation Details

for HPLC Application ID No.: 19714

Benzodiazepines from urine using Strata-X-Drug N on a Kinetex C18 2.6µm, 50x2.1mm

PRODUCT DESCRIPTION:

StrataTM-X-Drug N 100 µm Polymeric Reversed Phase, 100 mg / 6 mL, Tubes , 30/Pk

Order No.: 8B-S129-ECH

SOLID PHASE EXTRACTION (SPE) PROCEDURE:

Note: The solvent volumes shown below are for a 100 mg bed mass.

The solvent volumes will need to be adjusted for a smaller or larger bed mass.

Condition:

Load:

Wash:

Dry:

10 minutes @ full vacuum

Elute:

Final Prep and Analysis:

Evaporate eluant under a stream of N₂ gas at 50°C. Reconstitute with 1mL of 35% methanol

Inject: 0 µL on HPLC Mass Spectrometer (MS) @ amu (ambient)

ANALYTES:	Spiked Conc. (ng/mL)	Log P	pKa	% Rec	%RSC (n=0)
1 alpha-Hydroxyalprazolam	300			92.4	
2 Oxazepam-d5	300			100	
3 Oxazepam	300			105.3	
4 Alprazolam-d5	300			100	
5 Alprazolam	300			98.4	
6 Nordiazepam-d5	300			100	
7 Nordiazepam	300			103.7	
8 Lorazepam-d4	300			100	
9 Lorazepam	300			101.1	
10 Clonazepam-d4	300			100	
11 Clonazepam	300			99.6	
12 Temazepam-d5	300			100	
13 Temazepam	300			104.6	
14 Diazepam-d5	300			100	
15 Diazepam	300			101.1	

Note: This method is designed as a convenient starting point for further investigation and can be tailored to meet your extraction goals.
Call your local Phenomenex Representative for assistance in method development and optimization techniques.

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com.cn

info@phenomenex.com