

Bivalirudin on Synergi 4um Hydro-RP

Column: Synergi™ 4 µm Hydro-RP 80 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4375-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.1% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	40	40	60
	3	45	20	80
	4	55	80	20

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm ID, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

19728

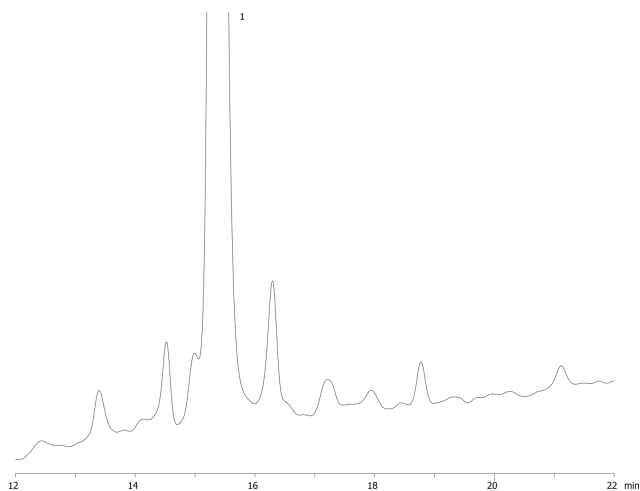

ANALYTES:

1 Bivalirudin

