

Intact Antibody on Aeris 3.6µm XB-C8 100x4.6mm

Column: Aeris™ 3.6 µm WIDEPORE XB-C8 200 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4481-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.1% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	1	80	20
	3	15	40	60

Flow Rate: 1 mL/min

Col. Temp.: 70 °C


Detection: UV-Vis Abs.-Diode Array (PDA) @ 280 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C8 for 4.6mm ID, 3/Pk Part No.: AJ0-8771

- Holder Part No.: AJ0-9000

19776

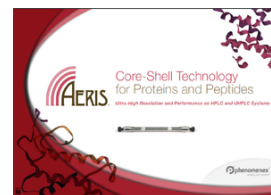

Products used in this application:

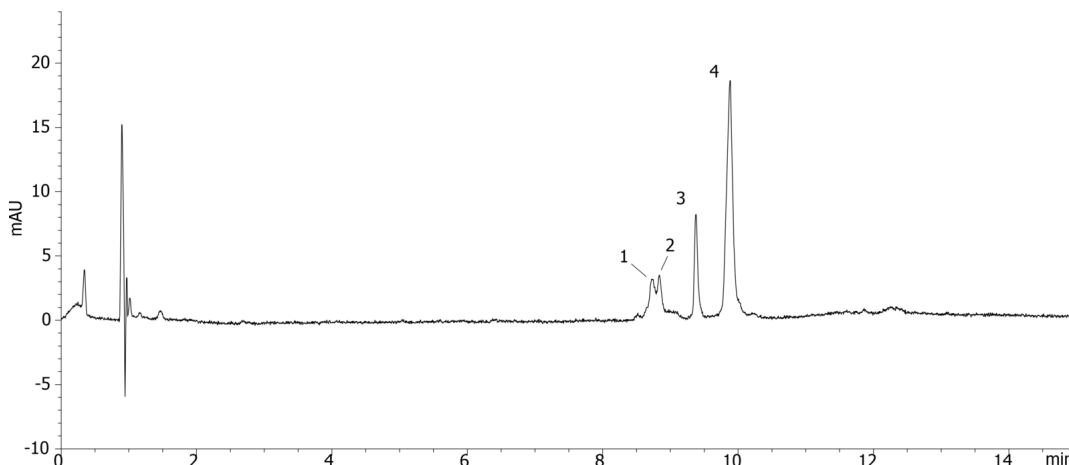

ANALYTES:

1 Mouse Monoclonal Antibody

