

Ethyl Sulfate (25ng/mL) on Synergi Hydro-RP (m/z 124.9 > 87.0)

Column: SynergiTM 2.5 µm Hydro-RP 100 Å, LC Column 100 x 3 mm, Ea

Dimensions: 100 x 3 mm ID

Order No: 00D-4387-Y0

Elution Type: Gradient

Eluent A: 10mM Ammonium formate

Eluent B: 50:50 Methanol: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	4	20	80
	3	4.1	100	0
	4	7	100	0

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

Detector Info: <a target="_blank"

19781

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20research&utm_source=phenomenex&utm_medium=referral">SCIEX<

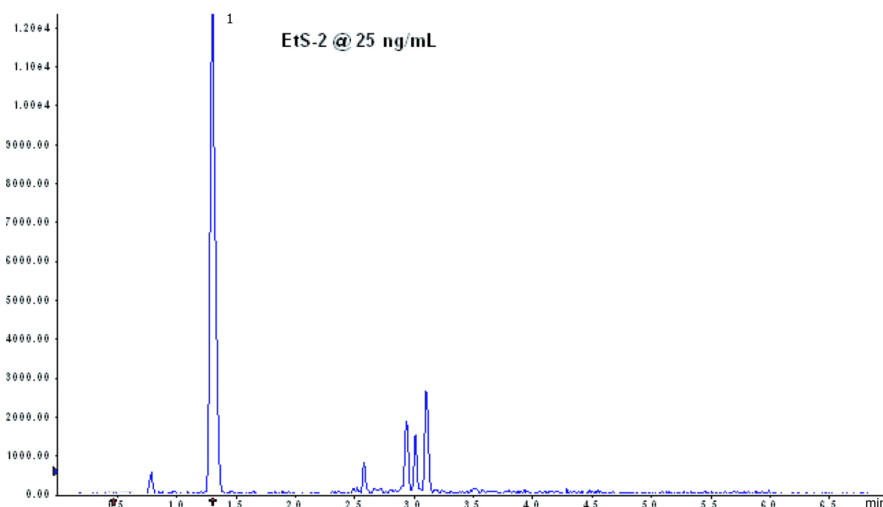

ANALYTES:

1 EtS (Ethyl sulfate)

