

Steviol Glycoside standards on Kinetex 2.6u XB-C18 150x4.6mm

Column: Kinetex® 2.6 µm XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4496-E0

Elution Type: Isocratic

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	68	32

Flow Rate: 1.5 mL/min

Col. Temp.: 70 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768


- Holder Part No.: AJ0-9000



Products used in this application:



19797

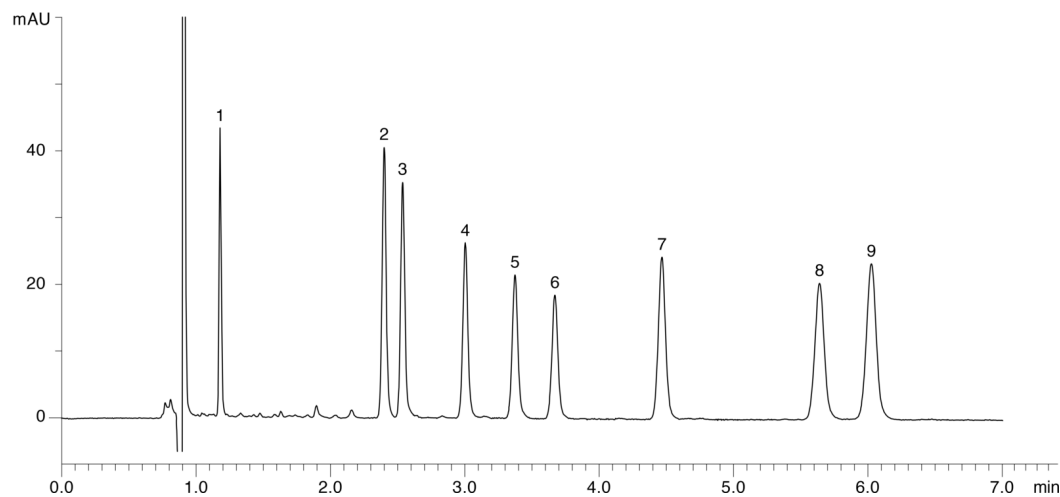

ANALYTES:

- 1 Rebaudioside D
- 2 Rebaudioside A
- 3 Stevioside
- 4 Rebaudioside F
- 5 Rebaudioside C
- 6 Dulcoside A
- 7 Rubusoside
- 8 Rebaudioside B
- 9 Steviolbioside

