

High Rebaudioside A sample on Kinetex 2.6u XB-C18 150x4.6mm**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4496-E0**Elution Type:** Isocratic**Eluent A:** Water with 0.1% TFA**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	68	32

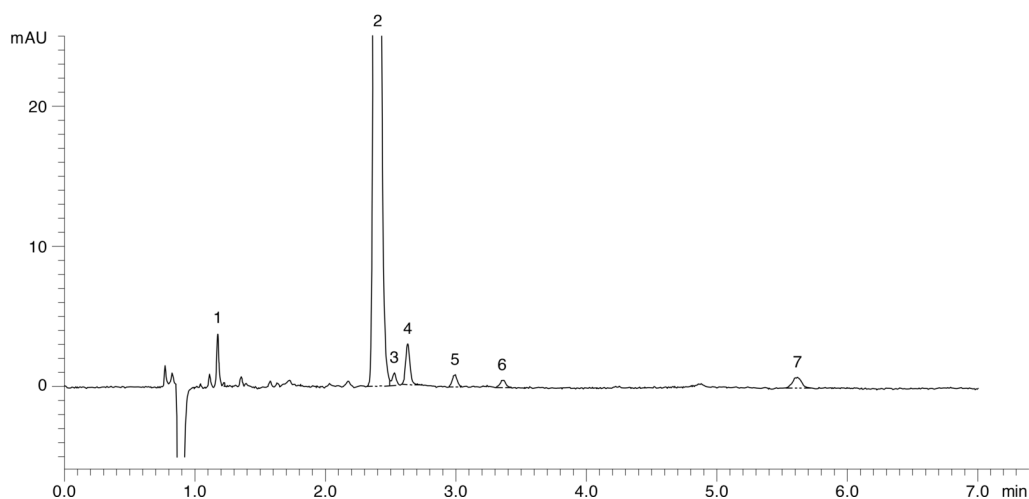
Flow Rate: 1.5 mL/min**Col. Temp.:** 70 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Analyst Note:** SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768

- Holder Part No.: AJ0-9000

**Products used in this application:**

19798

**ANALYTES:**

- 1 Rebaudioside D
- 2 Rebaudioside A
- 3 Stevioside
- 4 unknown
- 5 Rebaudioside F
- 6 Rebaudioside C
- 7 Rebaudioside B

