

GFC Standard on BioSep-sSEC-3000 (AL0-3042)

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Phosphate Buffer pH 6.8 with 0.025% NaN₃

Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, GFC 3000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4488

- Holder Part No.: KJ0-4282

19807

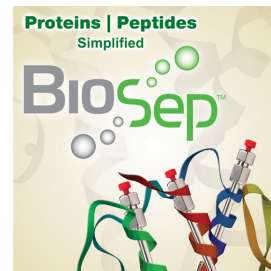

Products used in this application:

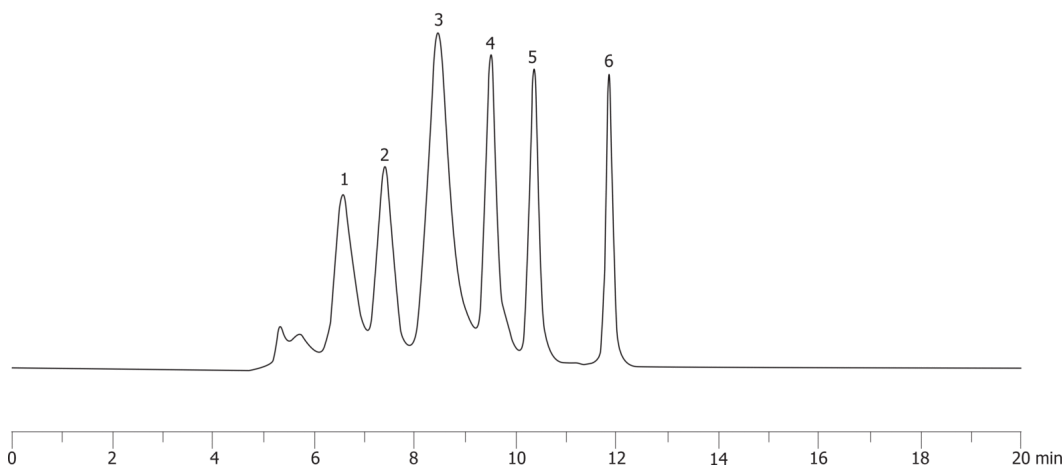

ANALYTES:

- 1** Bovine Thyroglobin
Retention Time: 6.56 min
- 2** IgA
Retention Time: 7.39 min
- 3** IgG
Retention Time: 8.45 min
- 4** Ovalbumin
Retention Time: 9.48 min
- 5** Myoglobin
Retention Time: 10.34 min
- 6** Uridine
Retention Time: 11.83 min

