

Aeris 3.6u C4 150 x 2.1mm ID QC Chromatogram (AL0-8931)

Column: Aeris™ 3.6 µm WIDEPORE C4 200 Å, LC Column 150 x 2.1 mm, Ea

Dimensions: 150 x 2.1 mm ID

Order No: 00F-4486-AN

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55

Flow Rate: 0.25 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C4 for 2.1mm ID, 3/Pk Part No.: AJ0-8899

- Holder Part No.: AJ0-9000

19810

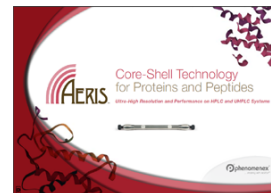

Products used in this application:

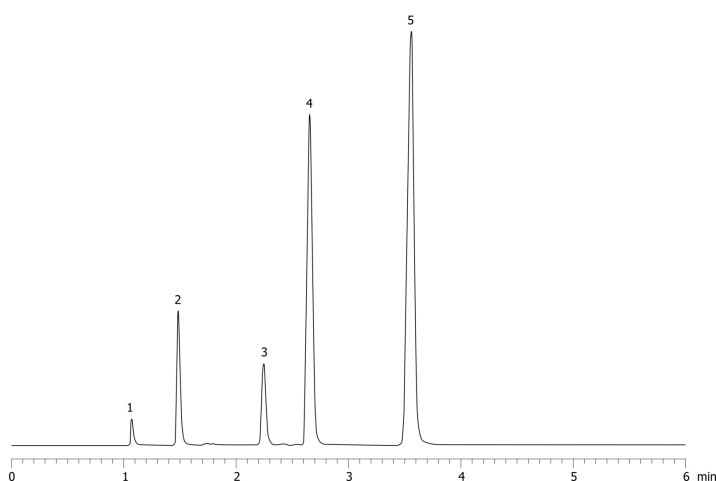

ANALYTES:

- 1** Uracil
Retention Time: 1.072 min
- 2** Acetophenone
Retention Time: 1.485 min
- 3** Toluene
Retention Time: 2.246 min
- 4** Naphthalene
Retention Time: 2.656 min
- 5** Acenaphthene
Retention Time: 3.558 min

