

Proprietary Monoclonal Antibody (Dyax) Aeris WIDEPORE 3.6u XB-C18 250x4.6mm

Column: Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4482-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile/IPA 50/50 with 0.1% TFA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	60	35	65

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 4.6mm ID, 3/Pk Part No.: AJ0-8769

- Holder Part No.: AJ0-9000

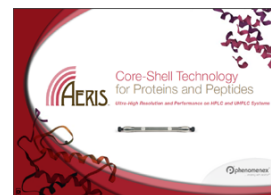

Products used in this application:

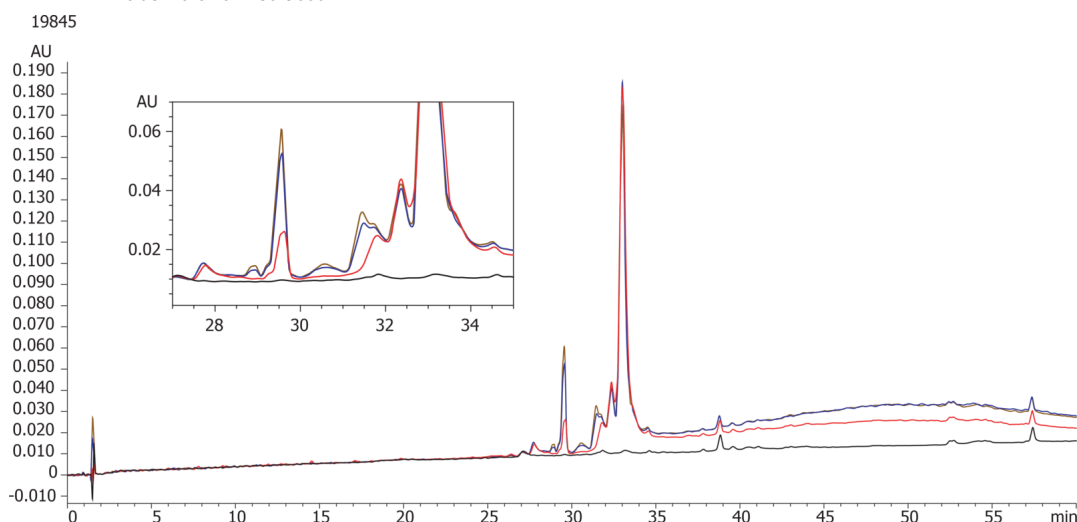

ANALYTES:

- 1 Monoclonal antibody / clipped variants

