

HPLC Application

ID No.: 19865

alpha-1-Acid Glycoprotein Trypsin Digest-Aeris PEPTIDE XB-C18 3.6µ 150x4.6mm

Column: Aeris™ 3.6 µm PEPTIDE XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4507-E0

Elution Type: Gradient

Eluent A: 0.1%TFA in Water


Eluent B: 0.1%TFA in ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	0	97	3
	3	3	97	3
	4	33	50	50
	5	35	5	95

Flow Rate: 1.2 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

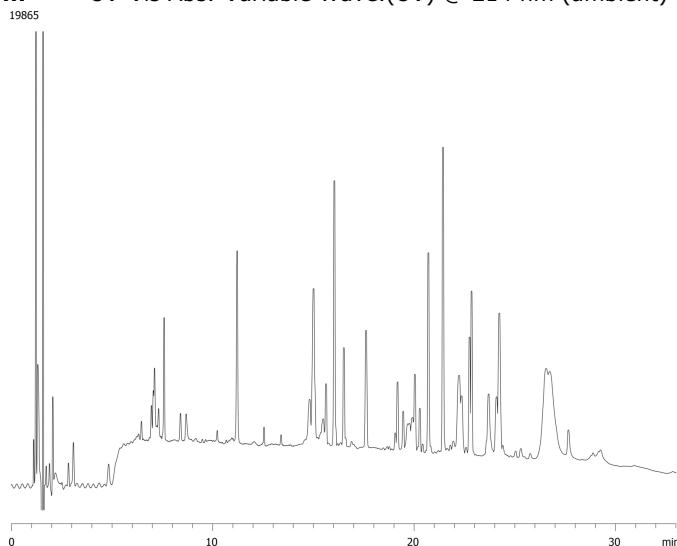

ANALYTES:

- 1 alpha-1-Acid Glycoprotein

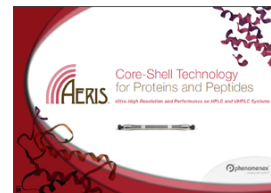

Products used in this application:

