

Protein Deamidation Aeris WIDEPORE 3.6um XB-C8 150x4.6mm

Column: Aeris™ 3.6 µm WIDEPORE XB-C8 200 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4481-E0

Elution Type: Gradient

Eluent A: 0.1% TFA/ Water

Eluent B: 0.085%TFA/ Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	3	10	35	65

Flow Rate: 1.2 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)



Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C8 for 4.6mm ID, 3/Pk Part No.: AJ0-8771

- Holder Part No.: AJ0-9000

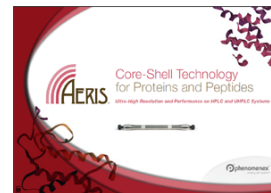

Products used in this application:

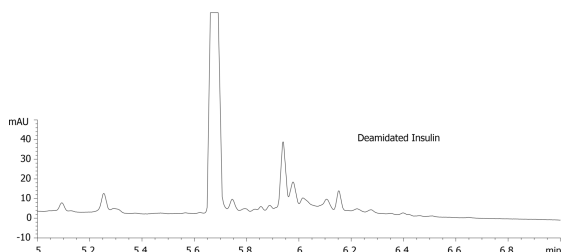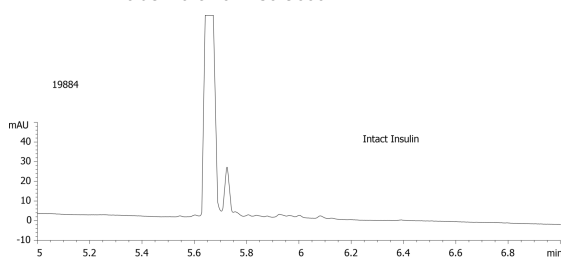

ANALYTES:

1 Intact Insulin

Retention Time: 5.6 min

