

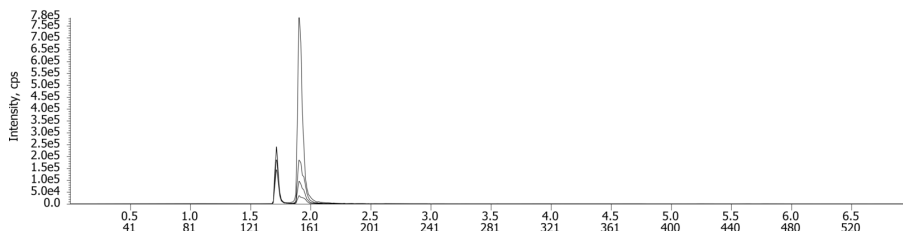
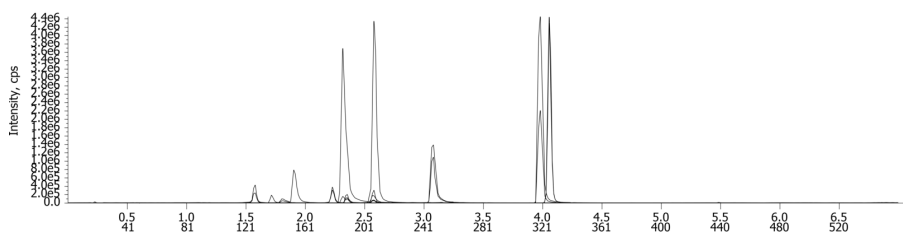
Opiates using Kinetex 2.6u XB-C18 50x2.1mm**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4496-AN**Elution Type:** Gradient**Eluent A:** 10mM ammonium formate**Eluent B:** 50/50 Methanol/Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** <a target="_blank"

19900

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**Products used in this application:****ANALYTES:****1** Opiates