

Security Guard ULTRA - No guard - RNAase A - on Aeris XB-C18 150x4.6mm

Column: Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4482-E0

Elution Type: Gradient

Eluent A: 0.1% TFA/ Water

Eluent B: 0.085% TFA/ Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	97	3
	2	3	97	3
	3	33	35	65

Flow Rate: 1.2 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 4.6mm ID, 3/Pk Part No.: AJ0-8769


- Holder Part No.: AJ0-9000

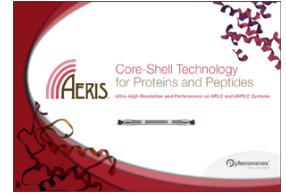

Products used in this application:



19923

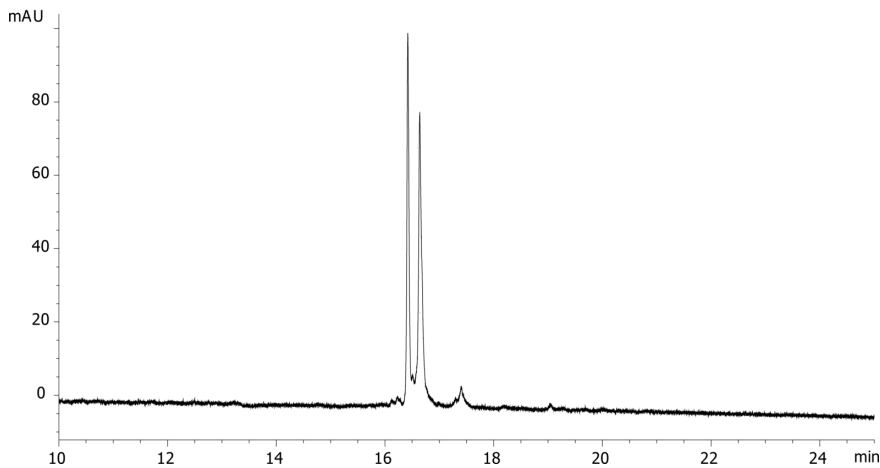

ANALYTES:

1 Intact RNase A

Retention Time: 16.42 min

2 Reduced RNase A

Retention Time: 16.639 min

