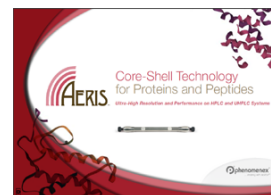
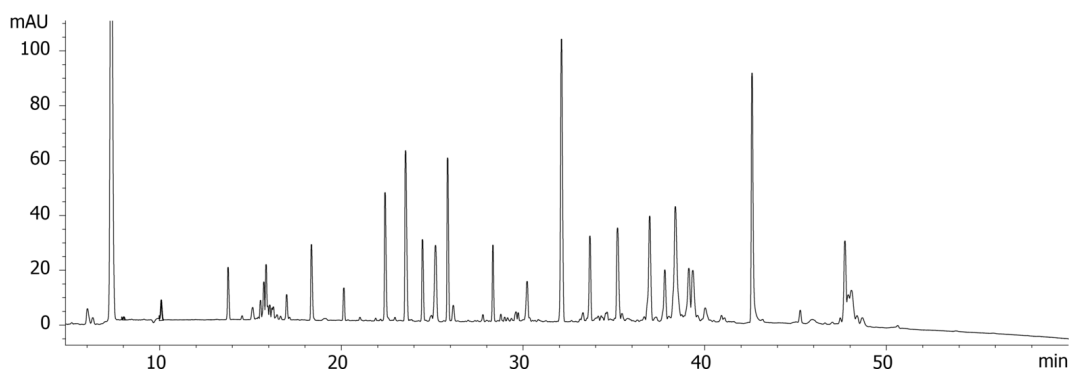


Lab Productivity Comparison Aeris PEPTIDE 3.6u XB-C18 150x4.6mm for Hu Ig-G**Column:** AerisTM 3.6 μ m PEPTIDE XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4507-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% TFA**Eluent B:** Acetonitrile with 0.1% TFA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	97	3
	2	5	97	3
	3	59	55	45

Flow Rate: 1.5 mL/min**Col. Temp.:** 40 °C**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)

19927

**Products used in this application:****ANALYTES:****1** Human IgG