

Phthalates in Foods by LC/MS /MS on Synergi Polar-RP 2.5µm, 100x2.0 mm

Column: Synergi™ 2.5 µm Polar-RP 100 Å, LC Column 100 x 2 mm, Ea

Dimensions: 100 x 2 mm ID

Order No: 00D-4371-B0

Elution Type: Gradient

Eluent A: 0.1% Formic acid in water

Eluent B: 0.1% Formic acid in 50:50 Acetonitrile: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	1	50	50
	3	6	50	50
	4	11	15	85
	5	14	15	85
	6	14.1	60	40
	7	19	60	40



Flow Rate: 0.4 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)


Detector Info: <a target="_blank"

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC Phenyl for 2.1mm ID Columns, 3/Pk Part No.: AJ0-8788

- Holder Part No.: AJ0-9000

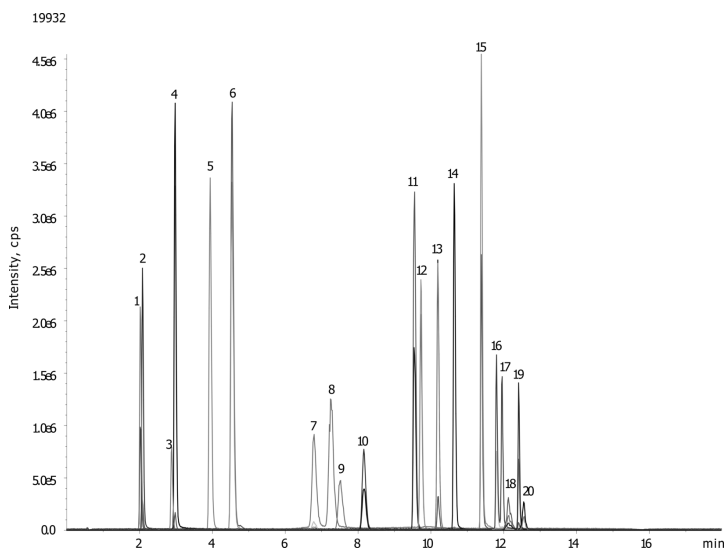

Phthalates in Foods by LC/MS /MS on Synergi Polar-RP 2.5µm, 100x2.0 mm

ANALYTES:

- 1** bis(2-Methoxyethyl)phthalate
- 2** Dimethyl phthalate
- 3** bis(2-Ethoxyethyl)phthalate
- 4** Diethyl phthalate
- 5** Diallyl phthalate
- 6** Di-n-propyl phthalate
- 7** Diisobutyl phthalate
- 8** Di-n-butyl phthalate
- 9** bis(2-n-Butoxyethyl)phthalate
- 10** Dibenzylbutyl phthalate
- 11** Dipentyl phthalate
- 12** Dicyclohexyl phthalate
- 13** Di-4-methyl-2-pentyl phthalate
- 14** Dihexyl phthalate
- 15** Diheptyl phthalate
- 16** bis(2-Ethylhexyl)phthalate
- 17** Di-n-octyl phthalate
- 18** Diisononyl phthalate
- 19** Di-n-nonyl phthalate
- 20** Diisodecyl phthalate

