

Pharmaceutical Drug Mix on Kinetex C18 2.6u 75 x 2.1mm @ 1.1 mL/min(Quotient Data)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 75 x 2.1 mm, Ea

Dimensions: 75 x 2.1 mm ID

Order No: 00C-4462-AN

Elution Type: Gradient

Eluent A: 0.1% formic acid in acetonitrile

Eluent B: 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	5	95
	2	0.25	5	95
	3	1.8	90	10
	4	1.9	90	10
	5	1.91	5	95
	6	2.5	5	95


Flow Rate: 1.1 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)



Products used in this application:

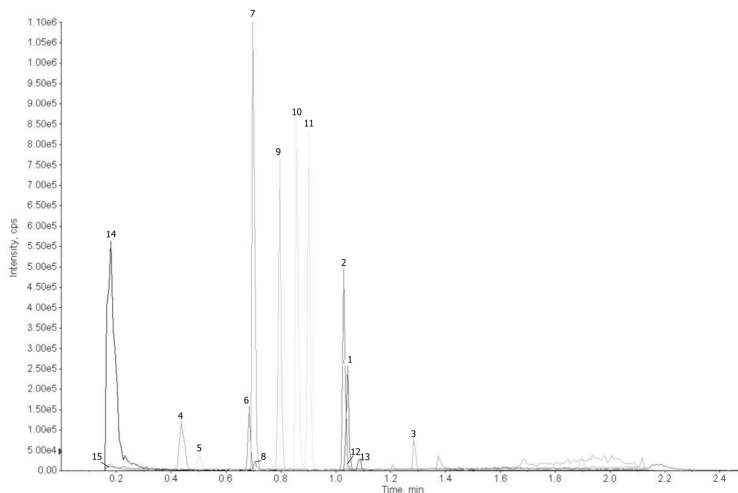

Pharmaceutical Drug Mix on Kinetex C18 2.6u 75 x 2.1mm @ 1.1 mL/min(Quotient Data)

ANALYTES:

- 1** Haloperidol
- 2** Diltiazem
- 3** Terfenadine
- 4** Cimetidine
- 5** Acetaminophen
- 6** Sulfathiazole
- 7** Pindolol
- 8** Quinidine
- 9** Acebutolol
- 10** Chlorpheniramine
- 11** Triprolidine
- 12** Prednisolone
- 13** Nortriptyline
- 14** 2-Hydroxy-5-methylbenzaldehyde
- 15** Hexanophenone

