

USP Norethindrone, Progesterone and Mestranol on Kinetex 2.6u C8 100x2.10mm on Shimadzu Nexera

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 100 x 2.1 mm, Ea

Dimensions: 100 x 2.1 mm ID

Order No: 00D-4497-AN

Elution Type: Isocratic

Eluent A: 50:50 Acetonitrile:Water

Gradient

Step No.	Time (min)	Pct A
1	10	100

Profile:


Flow Rate: 0.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)
19950



Products used in this application:

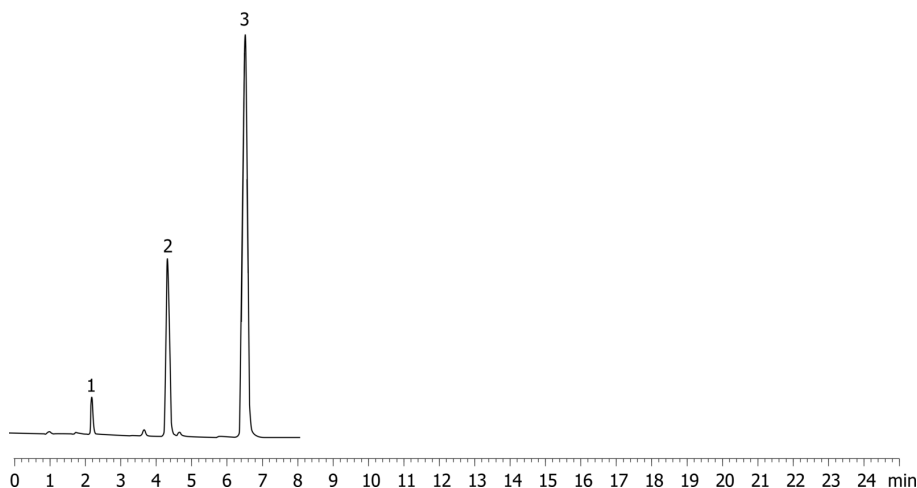

ANALYTES:

1 Norethindrone

Retention Time: 1.491 min

2 Progesterone

Retention Time: 2.968 min

3 Mestranol

Retention Time: 4.473 min

