

USP Norethindrone, Progesterone and Mestranol on Kinetex 2.6u C8 100x4.6mm on Shimadzu Nexera

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4497-E0

Elution Type: Isocratic

Eluent A: 50:50 Acetonitrile:Water

Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 25 100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)

19964



Products used in this application:

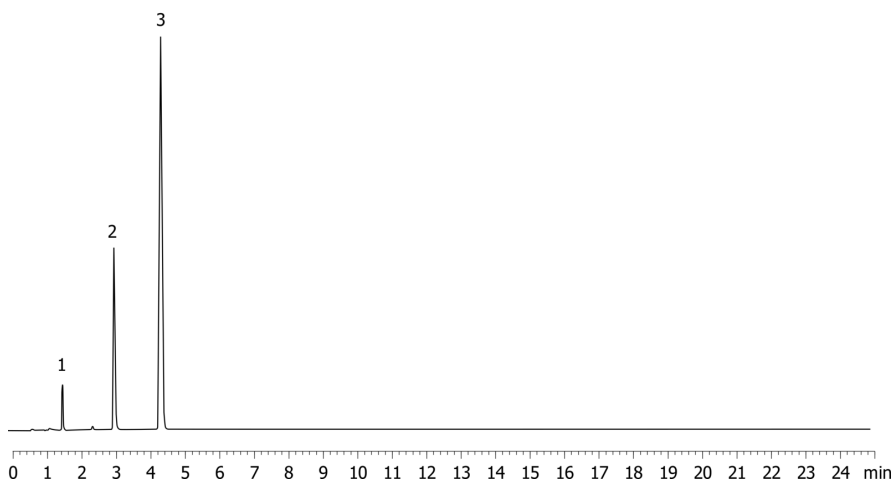

ANALYTES:

1 Norethindrone

Retention Time: 2.115 min

2 Progesterone

Retention Time: 4.326 min

3 Mestranol

Retention Time: 6.313 min

