

Column Guard Competitive Test: SecurityGuard ULTRA on Kinetex 1.7um XB-C18 50x2.1mm ID

Column: Kinetex® 1.7 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4498-AN

Elution Type: Isocratic

Eluent A: ACN

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50


Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)
19974



Products used in this application:

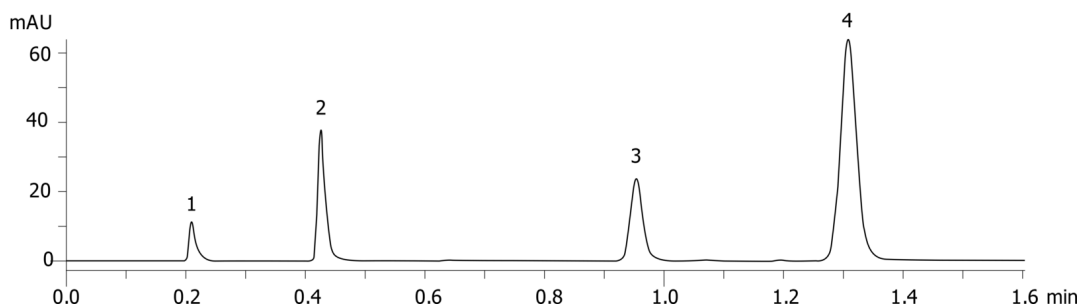

ANALYTES:

- 1 Acetophenone
- 2 Benzene
- 3 Toluene
- 4 Naphthalene

