

Isocratic Chromatogram of Kinetex with SecurityGuard ULTRA

Column: Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4496-AN

Elution Type: Isocratic

Eluent A: 50

Eluent B: 50

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0


Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)
19977



Products used in this application:

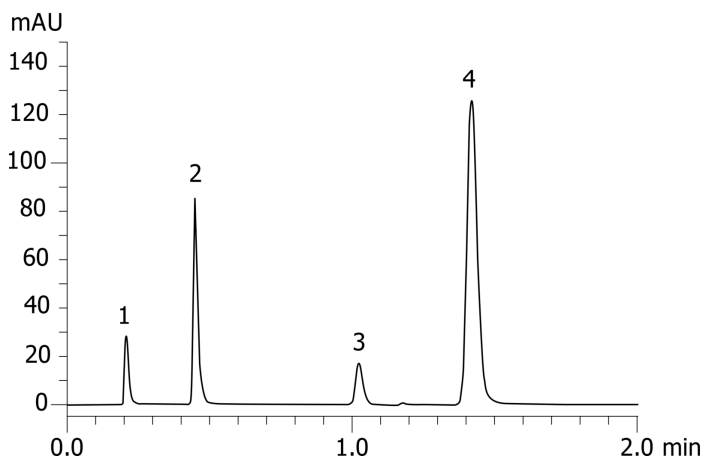

ANALYTES:

- 1** Uracil
Retention Time: 0.208 min
- 2** Acetophenone
Retention Time: 0.45 min
- 3** Toluene
Retention Time: 1.023 min
- 4** Naphthalene
Retention Time: 1.42 min

