

Kinetex 2.6 µm XB-C18 50 x 4.6 Prodigy Test With SecurityGuard ULTRA

Column: Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4496-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65


Flow Rate: 1.85 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ (ambient)
19980



Products used in this application:

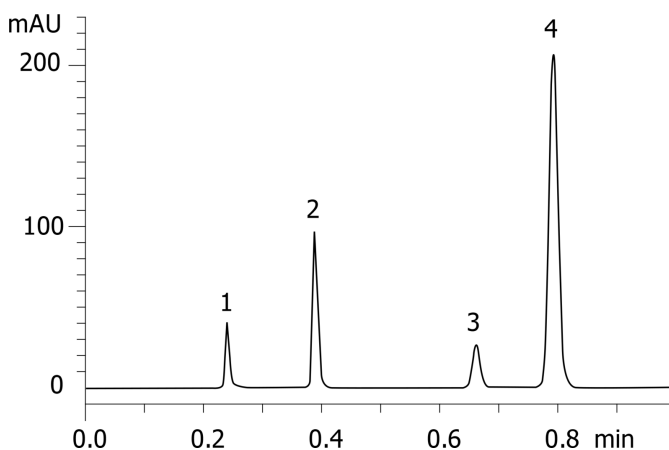

ANALYTES:

- 1** Uracil
Retention Time: 0.24 min
- 2** Acetophenone
Retention Time: 0.39 min
- 3** Toluene
Retention Time: 0.661 min
- 4** Naphthalene
Retention Time: 0.793 min

