

ABN+Aromatics KNX Phenyl-Hexyl 2.6μ 100Å

Column: Kinetex® 2.6 μm Phenyl-Hexyl 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4495-E0

Elution Type: Gradient

Eluent A: 0.1% Formic Acid/H₂O

Eluent B: 0.1% Formic Acid/MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	5	95

Flow Rate: 1.85 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19983



Products used in this application:

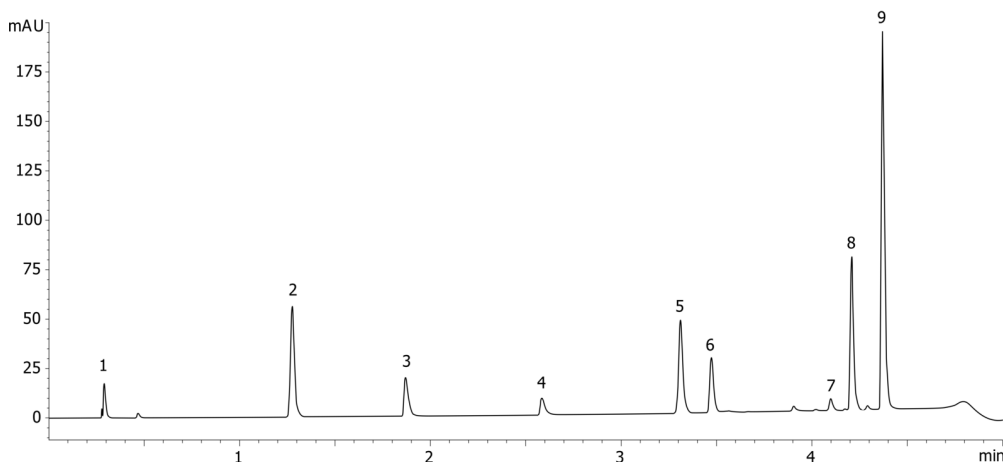

ANALYTES:

- 1** Uracil
Retention Time: 0.3 min
- 2** Acetaminophen
Retention Time: 1.28 min
- 3** Pindolol
Retention Time: 1.87 min
- 4** Chlorpheniramine
Retention Time: 2.59 min
- 5** 2-Hydroxy-5-methylbenzaldehyde
Retention Time: 3.31 min
- 6** 3-Methyl-4-nitrobenzaldehyde
Retention Time: 3.47 min
- 7** 1,5-Dinitronaphthalene
Retention Time: 4.1 min
- 8** 1,3-Dinitronaphthalene
Retention Time: 4.21 min
- 9** Anthracene
Retention Time: 4.37 min

