

Geometric Isomers XB-C18 2.6μ 100Å

Column: Kinetex® 2.6 μm XB-C18 100 Å, LC Column 50 x 4.6 mm, Ea
Dimensions: 50 x 4.6 mm ID
Order No: 00B-4496-E0
Elution Type: Isocratic
Eluent A: 20mM Potassium Phosphate pH 2.5
Eluent B: Acetonitrile
Gradient


Step No.	Time (min)	Pct A	Pct B
1	2	97	3

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

19988



Products used in this application:

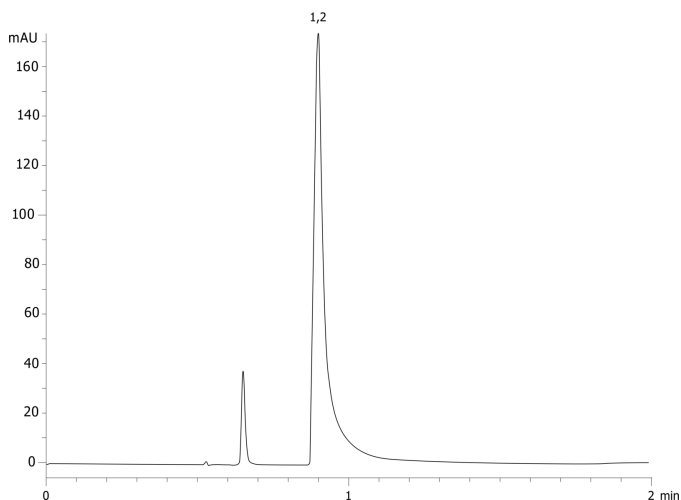

ANALYTES:

1 Fumaric acid / Maleic acid

