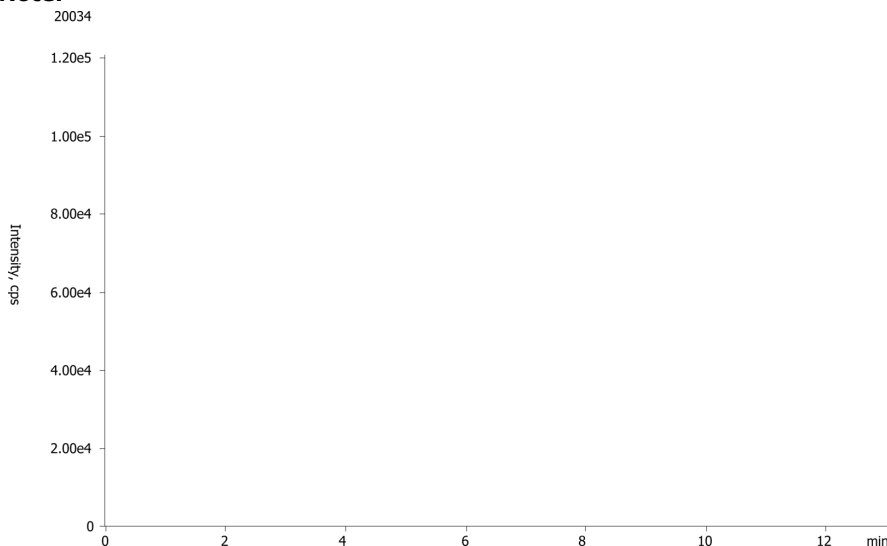


204 Pesticides in Food by iDQuant™ 10 Standards certified mixes**Column:** Synergi™ 2.5 µm Fusion-RP 100 Å, LC Column 50 x 2 mm, Ea**Dimensions:** 50 x 2 mm ID**Order No:** 00B-4423-B0**Elution Type:** Gradient**Eluent A:** Water/Methanol (90:10) + 5 mM ammonium formate**Eluent B:** Methanol/Water (90:10) + 5 mM ammonium formate

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	1	100	0
	3	15	0	100
	4	18	0	100
	5	18.05	100	0
	6	20	100	0

Flow Rate: 0.4 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX<,
Phenomenex SecurityGuard Guard Cartridge kit (P/N KJ0-4282) with C18 4 x 2.0 mm cartridge (P/N AJ0-4286)**ANALYTES:****1** junk