

Pharma ABN Test 15 Mix on Un-Optimized Agilent 1100 system using Kinetex 2.6u C18 100x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.2	95	5
	3	5	5	95
	4	6	5	95
	5	6.1	95	5

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.


- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768

- Holder Part No.: AJ0-9000

20049



Products used in this application:

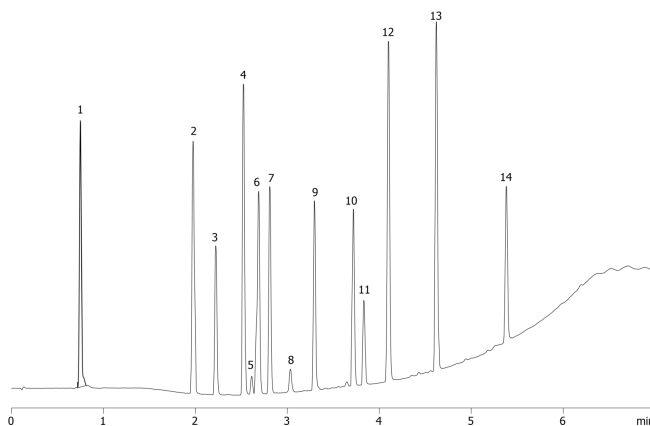

Pharma ABN Test 15 Mix on Un-Optimized Agilent 1100 system using Kinetex 2.6u C18 100x4.6mm

ANALYTES:

- 1** Pyridine
- 2** Acetaminophen
- 3** Sulfathiazole
- 4** Pindolol
- 5** Benzyl alcohol
- 6** Phenol
- 7** Acebutolol
- 8** Chlorpheniramine
- 9** Triprolidine
- 10** Prednisolone
- 11** 3-Methyl-4-nitrobenzoic acid
- 12** Nortriptyline
- 13** 2-Hydroxy-5-methylbenzaldehyde
- 14** Diflunisal
- 15** Hexanophenone

