

3,5-dinitro-N-[(1R)-1-phenylethyl]benzamide on Lux 3µm Cellulose-2 in NP

Column: Lux® 3 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4456-E0

Elution Type: Isocratic

Eluent A: Hexane / IPA (60:40) DEA (0.1%)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

20275



Products used in this application:

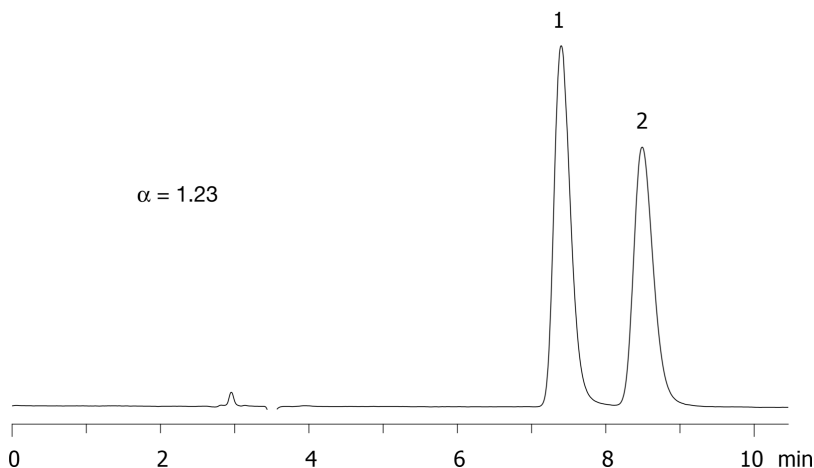

ANALYTES:

1 3,5-Dinitro-N-[(1R)-1-phenylethyl]benzamide

Retention Time: 7.4 min

2 3,5-Dinitro-N-[(1R)-1-phenylethyl]benzamide

Retention Time: 8.49 min

