

Suprofen on Lux 5µm Cellulose-3 in RP

Column: Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4493-E0

Elution Type: Isocratic

Eluent A: ACN / 0.1% Formic Acid (40:60) None

Gradient

Step No. 1 **Time (min)** 0 **Pct A** 100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622


- Holder Part No.: KJ0-4282



Products used in this application:



20294

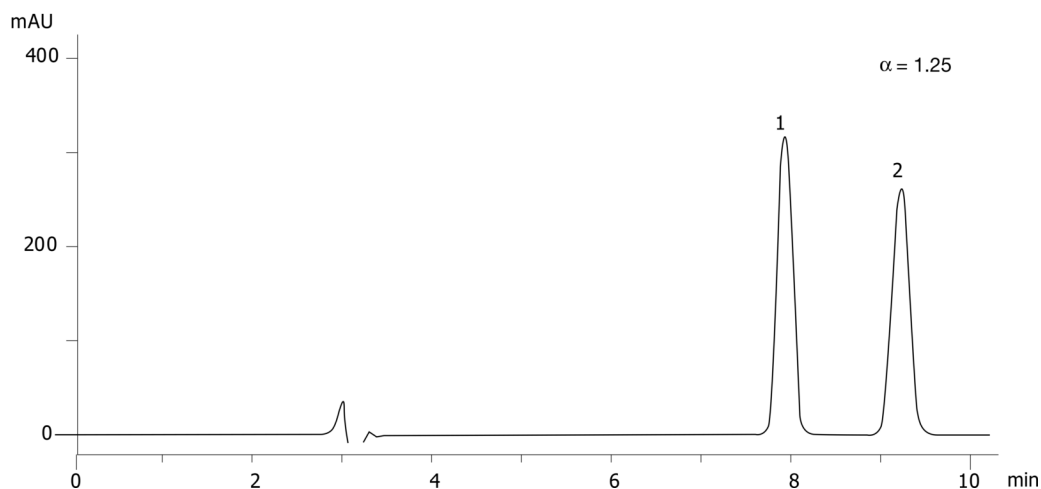

ANALYTES:

1 Suprofen

Retention Time: 7.93 min

2 Suprofen

Retention Time: 9.22 min

