

## Overlay: Yarra 3µm SEC-3000 and BioSep-SEC-S 3000 Low Salt PMQC Standard

**Column:** Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4513-K0

**Elution Type:** Isocratic


**Eluent A:** 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% Na<sub>3</sub>

| Gradient Profile: | Step No. | Time (min) | Pct A |
|-------------------|----------|------------|-------|
|                   | 1        | 0          | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

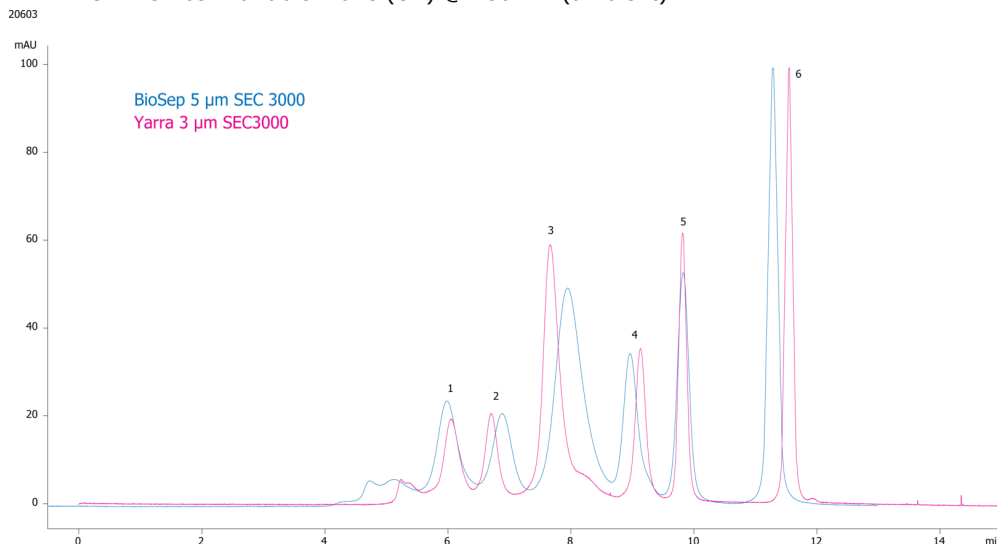

Products used in this application:



### ANALYTES:

- 1** Thyroglobulin  
Retention Time: 6 min
- 2** IgA  
Retention Time: 6.7 min
- 3** IgG  
Retention Time: 7.7 min
- 4** Ovalbumin  
Retention Time: 8.2 min
- 5** Myoglobin  
Retention Time: 9.8 min
- 6** Uridine  
Retention Time: 11 min

