


Yarra 3µm SEC-4000 PMQC Standard Low Salt Buffer

Column: Yarra™ 3 µm SEC-4000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4514-K0
Elution Type: Isocratic
Eluent A: 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% NaN₃
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20604

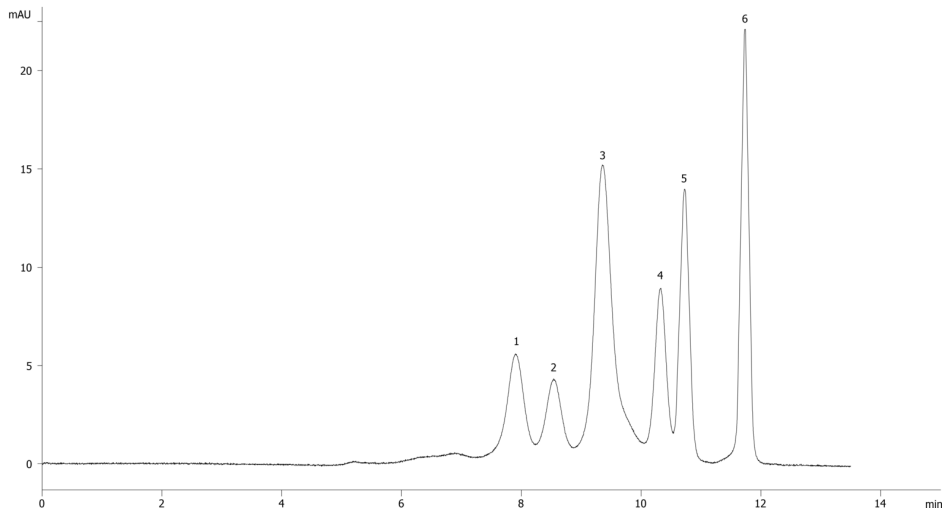

ANALYTES:

- 1** Thyroglobulin
Retention Time: 7.7 min
- 2** IgA
Retention Time: 8.5 min
- 3** IgG
Retention Time: 9.4 min
- 4** Ovalbumin
Retention Time: 10.3 min
- 5** Myoglobin
Retention Time: 10.7 min
- 6** Uridine
Retention Time: 11.7 min



Products used in this application:

