

BioSep-SEC-S 2000 High Salt Buffer MMQC Standard

Column: BioSep™ 5 µm SEC-s2000 145 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2145-K0

Elution Type: Isocratic

Eluent A: 50 mm Sodium Phosphate pH 6.8/ 0.3M Sodium Chloride

Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC-2000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4487

- Holder Part No.: KJ0-4282

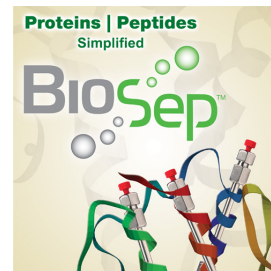

Products used in this application:

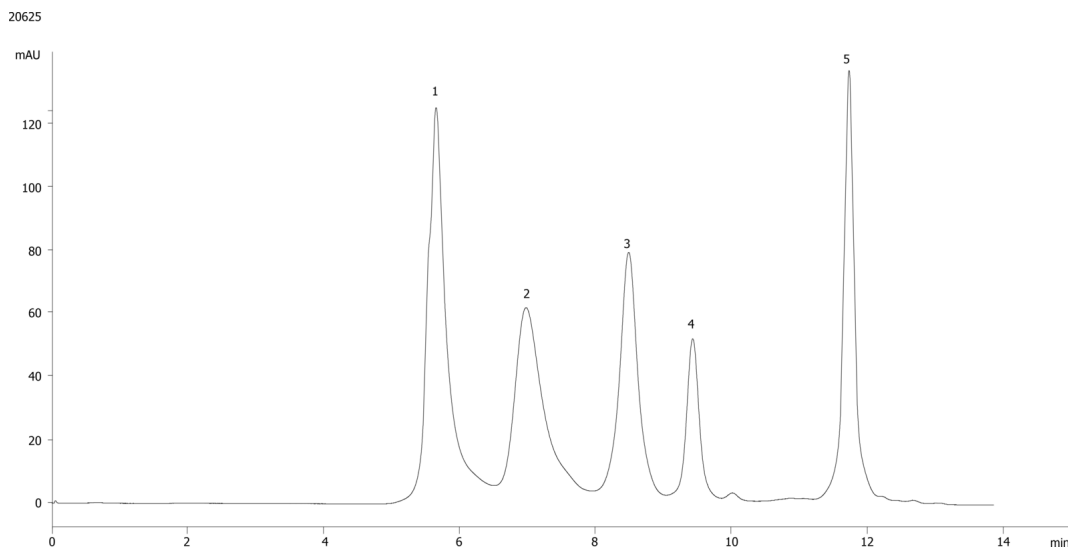

ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.7 min
- 2** Beta Amylase
Retention Time: 7 min
- 3** Ovalbumin
Retention Time: 8.5 min
- 4** Myoglobin
Retention Time: 9.5 min
- 5** Uridine
Retention Time: 11.7 min
- 6** Uridine
Retention Time: 11.3 min

