


Yarra 3µm SEC-2000 high salt buffer MMQC Standard

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4512-K0
Elution Type: Isocratic
Eluent A: 50mM Sodium Phosphate pH 6.8/ 0.3M Sodium Chloride
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20626

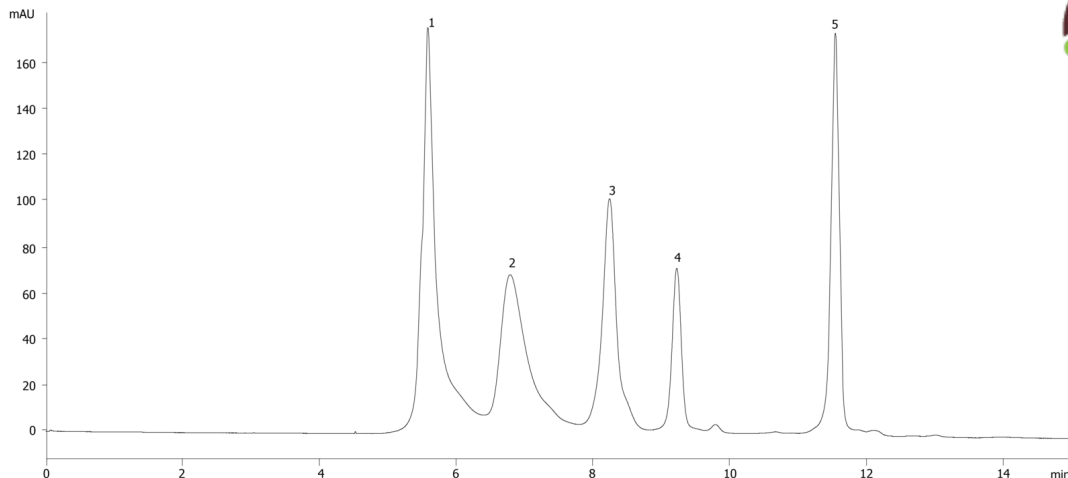

ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.7 min
- 2** Beta Amylase
Retention Time: 6.8 min
- 3** Ovalbumin
Retention Time: 8.3 min
- 4** Myoglobin
Retention Time: 9.2 min
- 5** Uridine
Retention Time: 11.5 min



Products used in this application:

