


Yarra 3µm SEC-3000 High Salt Buffer MMQC Standard

Column: Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4513-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate pH 6.8/ 0.3M Sodium Chloride
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20628

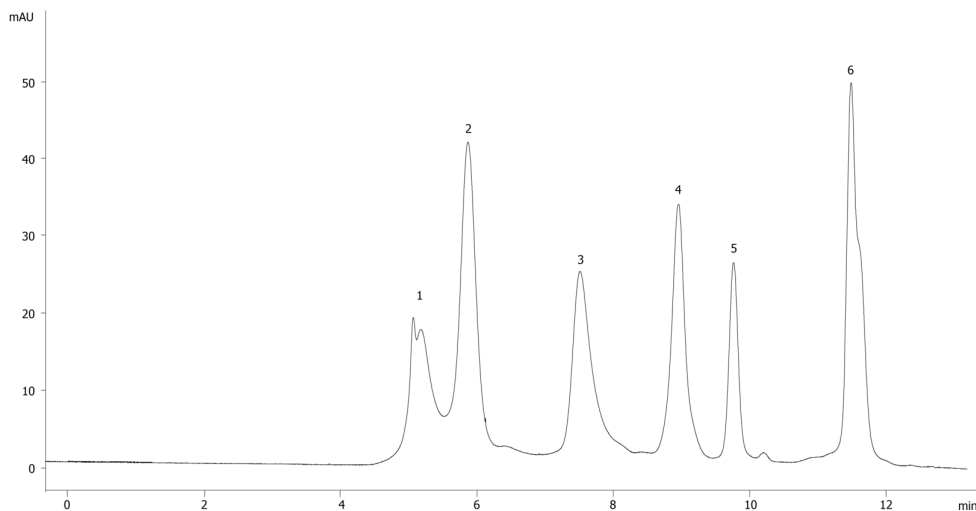

ANALYTES:

- 1** IgM
Retention Time: 5.5 min
- 2** Thyroglobulin
Retention Time: 6.3 min
- 3** Beta Amylase
Retention Time: 7.9 min
- 4** Ovalbumin
Retention Time: 9.3 min
- 5** Myoglobin
Retention Time: 10.1 min
- 6** Uridine
Retention Time: 11.7 min



Products used in this application:

