


Overlay: Yarra 3µm SEC-2000 versus TSKgel G2000 SWxl high salt buffer MMQC Standard

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4512-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride
Gradient Profile:

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20630

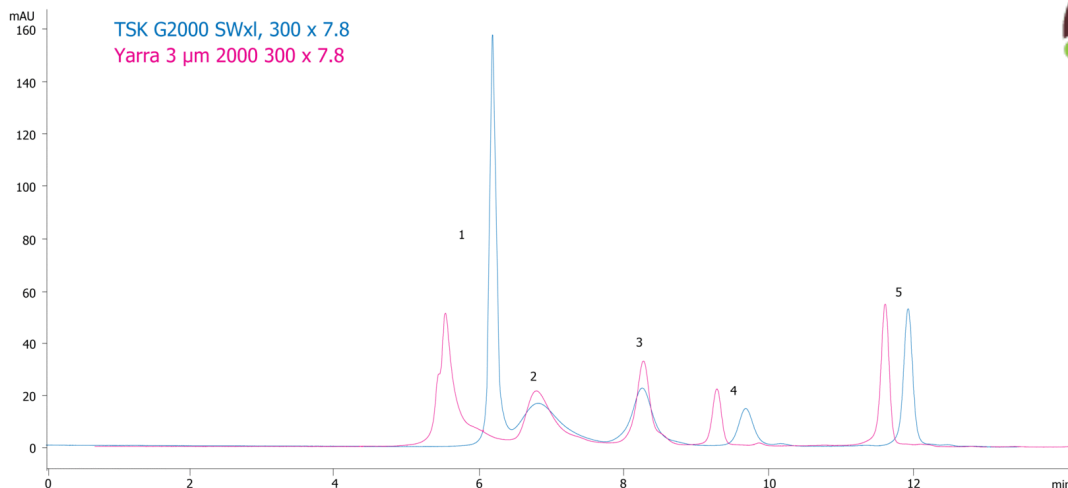

Products used in this application:



ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.7 min
- 2** Beta Amylase
Retention Time: 7 min
- 3** Ovalbumin
Retention Time: 8.5 min
- 4** Myoglobin
Retention Time: 9.5 min
- 5** Uridine
Retention Time: 11.7 min

