

Calibration curves for Yarra phases SEC-2000, SEC-3000, SEC-4000

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-4512-K0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Chloride pH 6.8

Gradient **Step No.** **Time (min)** **Pct A**
Profile: **1** 0 100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)
20638



Products used in this application:

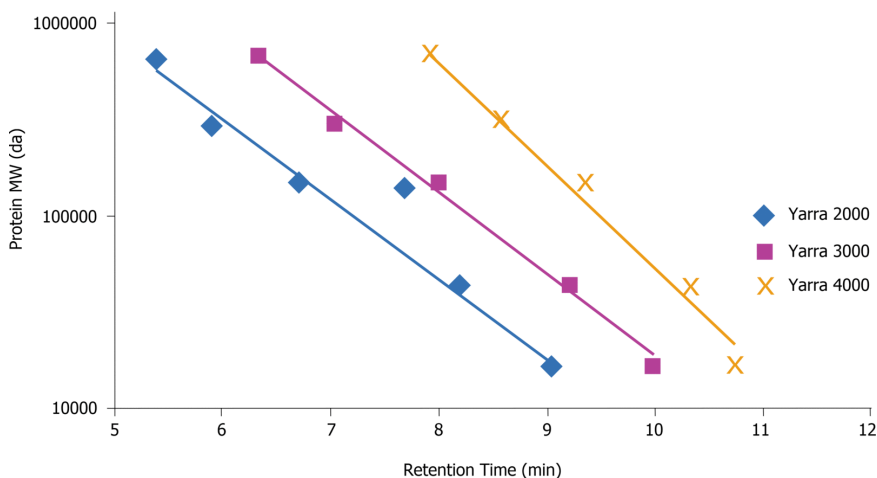

ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine

