


Gabapentin by LCMS using Kinetex 2.6u C18 50x2.1mm**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4462-AN**Elution Type:** Gradient**Eluent A:** 5mM ammonium formate**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1	97	3
	3	3	10	90
	4	3.1	97	3
	5	5	97	3

Flow Rate: 0.4 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX**Products used in this application:****ANALYTES:****1** Gabapentin

Retention Time: 1.44 min

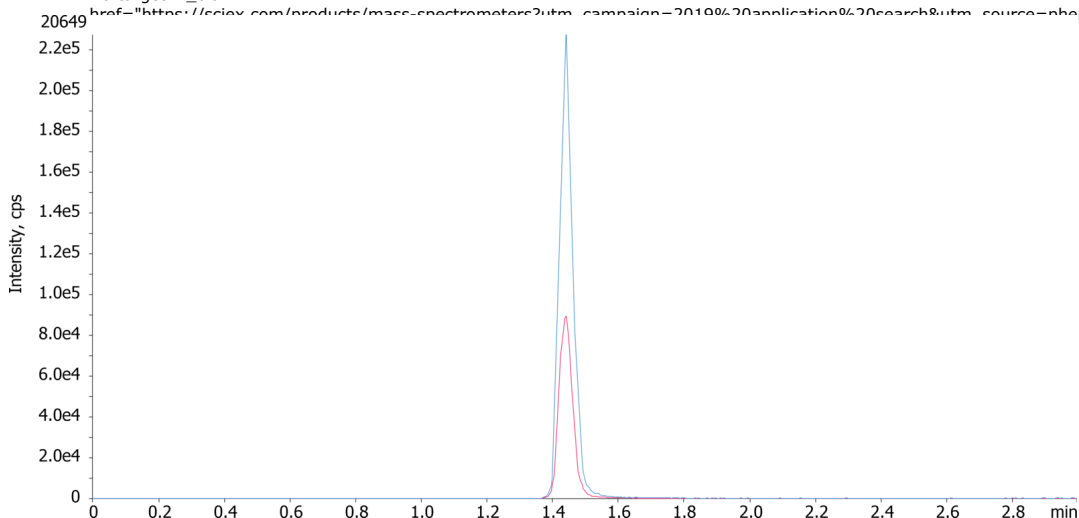