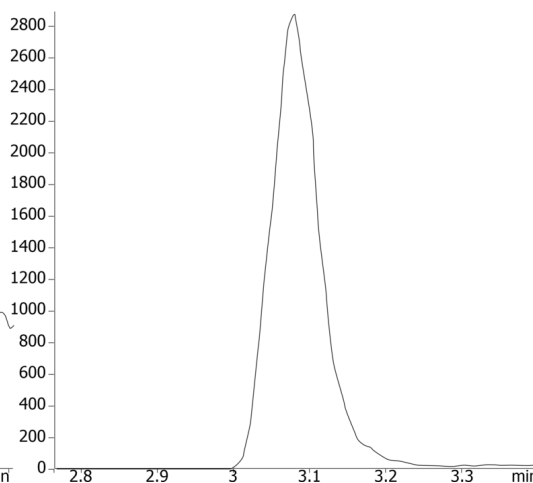
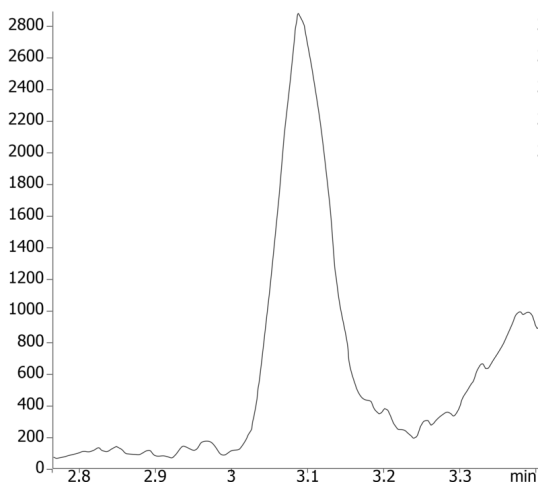


Estriol by LCMS using Gemini 3um C6-Phenyl 50x2mm**Column:** Gemini® 3 µm C6-Phenyl 110 Å, LC Column 50 x 2 mm, Ea**Dimensions:** 50 x 2 mm ID**Order No:** 00B-4443-B0**Elution Type:** Gradient**Eluent A:** H₂O + 0.5 mM ammonium hydroxide**Eluent B:** MeOH + 0.5 mM ammonium hydroxide

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	0.5	80	20
	3	3.5	10	90
	4	5	10	90
	5	5.1	80	20
	6	6	80	20

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX<

20654

**ANALYTES:****1** Estriol

Retention Time: 3.08 min

2 Estriol-d3

Retention Time: 3.08 min



Products used in this application:

