

Degraded Myoglobin on Jupiter 5u C18 compare to aeris

Column: Jupiter® 5 µm C18 300 Å, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4053-B0

Elution Type: Gradient

Eluent A: 0.1% TFA

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	50	50	50

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000



Products used in this application:

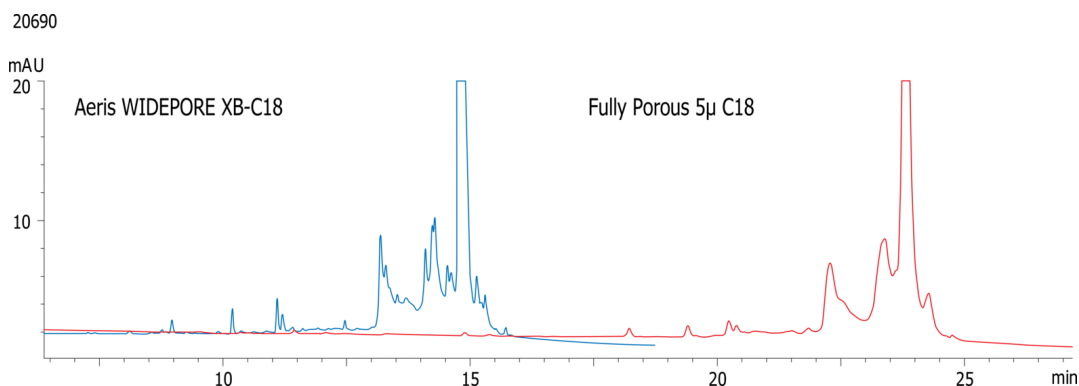

ANALYTES:

1 Myoglobin

