

Omeprazole on Lux 5 μ Cellulose-1 in NP Elution Mode (0.1% FA)

Column: Lux[®] 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: 0.1% FA IN HEXANE

Eluent B: 0.1% FA IN IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	20	40	60


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)
20701



Products used in this application:

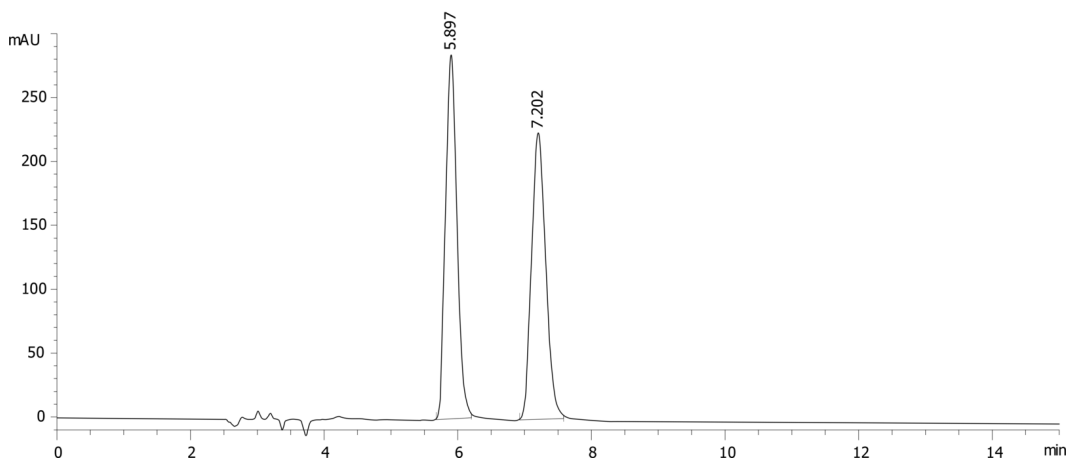

ANALYTES:

- 1 Omeprazole
- 2 Omeprazole

